YILKAL M E D I C A L





Nebulizer set

It makes the medicine in liquid form to an aerosol and provides the medicine to reach the respiratory tract during the patient's breathing with the help of the mask.

Features

Masks are designed in four different sizes according to the needs of newborn, infant, child and adult patients. Soft and compatible with the facial features masks with the nose clip covers the mouth and nose. Fixation elastic that surrounds the head. The oxygen connection hose is designed with multiple channels and is produced resistant to breakage and bending, so that the oxygen flow continues without obstruction. It is compatible with central and other gas supply systems with the standard connection connector.

Medicine Reservoir With the low residual rate calibrated reservoir, the nebulizer medicine loss is minimized. It can be opened and closed with a quarter turn with the the lid lock system that allows quick mounting to the reservoir. In ongoing treatments, medication can be added by opening the small cover on the top while the set is on the patient face. The reservoir delivers high performance aerosol providing optimum particle size at angles up to 90 degrees. Sealing lid to prevent spillage and 6cc capacity reservoir.

Quality

- * ISO 13485:2016
- * CE certificated

Shelf Life

5 Years



REF NO	PRODUCT NAME	TUBE	PACKING
102 020 1	Nebulizer Set Newborn	190 cm.	200 Pcs/Crt 13,20 kg. 37x58x38 cm
102 020 2	Nebulizer Set Infant	190 cm.	175 Pcs/Crt 12,55 kg. 37x58x38 cm
102 020 3	Nebulizer Set Pediatric	190 cm.	150 Pcs/Crt 11,41kg. 37x58x38 cm
102 020 4	Nebulizer Set Adult	190 cm.	125 Pcs/Crt 10,12 kg. 37x58x38 cm

B.02/24.05.21/REV0





Nebulizer Set

It makes the medicine in liquid form to an aerosol and provides the medicine to reach the respiratory tract during the patient's breathing with the help of the mouthpiece.

Designed in one size according to the needs of newborn, infant, pediatric and adult patients. Suitable for mouth anatomy. The oxygen connection hose is designed with multiple channels and is produced resistant to breakage and bending, so that the oxygen flow continues without obstruction. It is compatible with central and other gas supply systems with the standard connection

Medicine Reservoir With the low residual rate calibrated reservoir, the nebulizer medicine loss is minimized. It can be opened and closed with a quarter turn with the lid lock system that allows quick mounting to the reservoir. In ongoing treatments, medication can be added by opening the small cover on the top while the set is on the patient face. The reservoir delivers high performance aerosol providing optimum particle size at angles up to 90 degrees. Sealing lid to prevent spillage and 6cc capacity reservoir.

- Quality
 * ISO 13485:2016
- * CE Certificated

Shelf Life

5 Years



REF NO	PRODUCT NAME	TUBE	PACKING
103 020 1	Nebulizer Set With Mouthpiece	190 cm.	175 Pcs/Crt 10.36 kg.37x58x38 cm

B.01/24.05.21/REV0



Nebulizer Set

It aerosolizes the liquid medicine and allows it to reach the respiratory tract from the respiratory circuit through the T-piece.


Features

It is designed in universal sizes for use with ventilator circuits as a uniform. The oxygen connection hose is designed with multiple channels and is produced resistant to breakage and bending, so that the oxygen flow continues without obstruction. It is compatible with central and other gas supply systems with the standard connection connector.

Medicine Reservoir With the low residual rate calibrated reservoir, the nebulizer medicine loss is minimized. It can be opened and closed with a quarter turn with the lid lock system that allows quick mounting to the reservoir. In ongoing treatments, medication can be added by opening the small cover on the top while the set is on the patient face. The reservoir delivers high performance aerosol providing optimum particle size at angles up to 90 degrees. Sealing lid to prevent spillage and 6cc capacity reservoir.

- **Quality*** ISO 13485:2016
- * CE Certificated

Shelf Life

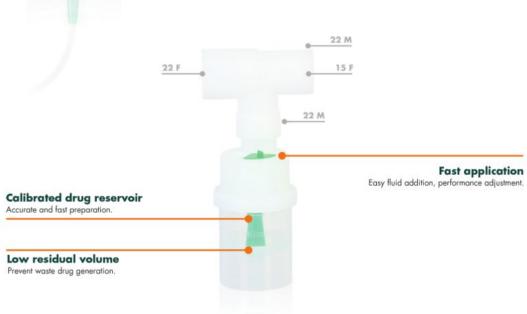

REF NO	PRODUCT NAME	TUBE	PACKING	
103 020 1	Nebulizer Set With T Piece	190 cm.	175 Pcs/Crt 11 32kg 37x58x38cm	





Nasal Oxygen Cannula
Provides easy and safe nasal oxygen support to patients with respiratory distress and decreased saturation.

It is designed for four sizes as newborn, infant, pediatric and adult according to the nose sizes and patient needs. Nosepieces have two open ends. These pieces and the oxygen hose do not give the patient a feeling of weight with their soft and flexible structure. All oxygen hose series are designed with multi-channels and produced resistant to breakage and bending so that the oxygen flow continues without obstruction. The cannula can be fixed in the desired position with the ring at the back of the head. It is compatible with central and other gas supply systems with the standard connection connector.

- **Quality*** ISO 13485:2016
- * CE certificated

Shelf Life

5 Years

Convenient use

Soft and flexible for patient comfort.

Fast and safe use

Non-slip easy fixing, increased patient comfort.

Ī	REF NO	PRODUCT NAME	TUBE	PACKING
ı	104 020 1	Nasal Oxygen Cannula Newborn	200 cm.	325 Pcs/Crt 11,00 kg. 37x58x38 cm
	104 020 2	Nasal Oxygen Cannula Infant	200 cm.	325 Pcs/Crt 11,20 kg. 37x58x38 cm
	104 020 3	Nasal Oxygen Cannula Pediatric	200 cm.	325 Pcs/Crt 11,30 kg. 37x58x38 cm
	104 020 4	Nasal Oxygen Cannula Adult	200 cm.	325 Pcs/Crt 11,40 kg. 37x58x38 cm



HIGH CONCENTRATION OXYGEN MASK



High Concentration Oxygen Mask

Provides oxygen support at high concentrations such as 80-99% in an easy and safe way, to patients with respiratory distress and decreased saturation rate.


Features

Masks are designed in two different sizes according to the needs of child and adult patients. Soft and compatible with the facial features masks with the nose clip covers the mouth and nose. Fixation elastic that surrounds the head. It has low resistance check valves to prevent rebreathing and allow exhaled gas to escape. Contains 750 ml oxygen reservoir. The oxygen connection hose is designed with multiple channels and is produced resistant to breakage and bending, so that the oxygen flow continues without obstruction. It is compatible with central and other gas supply systems with the standard connection connector.

- Quality * ISO 13485:2016
- * CE certificated

Shelf Life

5 Years

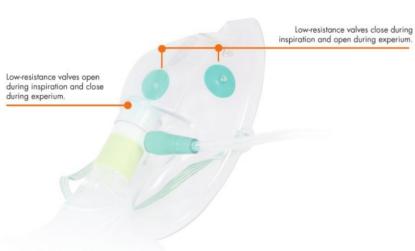

REF NO	PRODUCT NAME	TUBE	PACKING
112 020 3	High Concentration Oxygen Mask Pediatric	190 cm.	150 Pcs/Crt 9,90 kg. 37x58x38 cm
112 020 4	High Concentration Oxygen Mask Adult	190 cm.	125 Pcs/Crt 11,20 kg. 37x58x38 cm

B.06/24.05.21/REV0





Oxygen Mask

Provides oxygen support to patients with respiratory distress and decreased saturation rate.

Masks are designed in four different sizes according to the needs of newborn, infant, child and adult patients. Soft and compatible with the facial features masks with the nose clip covers the mouth and nose. Fixation elastic that surrounds the head. The oxygen connection hose is designed with multiple channels and is produced resistant to breakage and bending, so that the oxygen flow continues without obstruction. It is compatible with central and other gas supply systems with the standard connection connector.

- **Quality*** ISO 13485:2016
- * CE certificated

Shelf Life





Vaginal Speculum


It helps the gynaecologist to open the vagina. So that she can perform a vaginal examination and see the cervix comfortably. It also helps to obtain the necessary cervical cells for pop

Features

It has been designed for three sizes according to patient types as small, medium and large. It has a transparent design that provides easy use and offers high visual quality. 4 cm viewing aperture, fast locking at any desired angle, flexible lock system produced in different colors for all sizes. Strong design based on 40 N opening force. Rounded ergonomic tip that prevents tissue damage. Non-slip handle design that provides grip security.

- Quality * ISO 13485:2016
- * CE certificated

Shelf Life

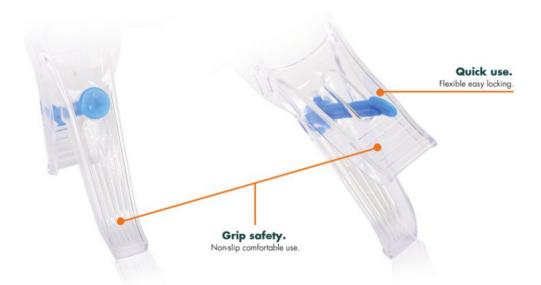

REF NO	PRODUCT NAME	LENGTH	PACKING
106 008 1	Vaginal Speculum Small	8-9 cm.	150 Pcs/Crt 7,16 kg. 37x58x38 cm
106 009 2	Vaginal Speculum Medium	9-10 cm.	140 Pcs/Crt 7,30 kg. 37x58x38 cm
106 010 3	Vaginal Speculum Large	10-11 cm.	125 Pcs/Crt 6,50 kg. 37x58x38 cm





Vaginal Speculum
It helps the gynaecologist to open the vagina. So that she can perform a vaginal examination and see the cervix comfortably. It also helps to obtain the necessary cervical cells for pop

It has been designed for two sizes according to patient types as medium and large. 6*3.5 cm viewing aperture, locking system on the handle that allows one-handed locking at the desired angle. Strong design based on 30 N opening force. Rounded ergonomic tip that prevents tissue damage. Non-slip handle design that provides grip security.

- Quality * ISO 13485:2016
- * CE certificated

Shelf Life



	REF NO	PRODUCT NAME	LENGTH	PACKING
Ī	106 009 5	Vaginal Speculum Medium	9-10 cm.	140 Pcs/Crt 7,23 kg 37x58x38 cm
	106 010 6	Vaginal Speculum Large	10-11 cm.	125 Pcs/Crt 7,06 kg 37x58x38 cm